

BS-2011

Bone Saw

(Rotate Mini)



Application and Feature

Electric surgical drill & saw is one of power tools widely used in orthopedic operation. Gradually replaced traditional industrial drill and saw, it can be autoclaved up to 135 °C high temperature and 1.6Mpa high pressure. Some of bacteria can be killed by high temperature and pressure, fumigating cannot do. Meanwhile, it can shorten sterilization time. Under high temperature and high pressure, just need 20 minutes, win the valuable time for the patient. The device is optional, good looking, small size, light weight, no wires, easy to carry, safe operation.

Technical Data

Item	Content	Parameter
Handpiece	Working voltage	7.2V
	Output power	≥50W
	Unit rise in temperature	≤20°C
Charger	Input voltage	AC220V,50Hz
	Output power	7.2V
	Input power	30VA
Battery	Voltage	7.2V
	Charging time	About 3 hours
Saw	autoclavable	up to 135 °C to do sterilization
	Frequency	0-18000 rpm
	Non-load Noise	≤45dB

Composition

Battery operated surgical drill & saw is made of main hand piece, charger, battery. The whole outer skin of the hand piece is sealed. Charger is not a medical electric device.

Standard Configuration

One saw unit
One charger
Two rechargeable batteries
Two aseptic battery transfer rings
Two saw blades
One instruction manual
One Aluminum case

Accessories

Accessories



Battery of 7.2V



Charger




Aseptic Battery Transfer Rings



Full Set



Sterilization Case

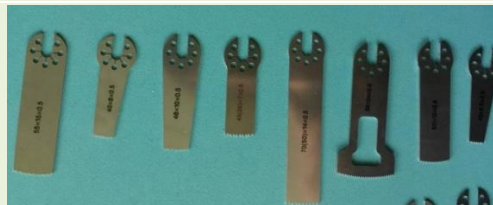

Saw blades